

Contents

1. Description
 - 1.1 Background information
 - 1.2 Applications
 - 1.3 Reagent and instrument requirements
2. Protocol
 - 2.1 Immunofluorescent staining
 - 2.2 Calibration of the MACSQuant Instrument
 - 2.3 Compensation
3. Example of a compensation using the MACS Comp Bead Kit

Warnings

Reagents contain sodium azide. Under acidic conditions sodium azide yields hydrazoic acid, which is extremely toxic. Azide compounds should be diluted with running water before discarding. These precautions are recommended to avoid deposits in plumbing where explosive conditions may develop.

1. Description

This product is for research use only.

Components	5 mL MACS Comp Beads – anti-REA 5 mL MACS Comp Beads – blank
Capacity	For 100 tests.
Product format	MACS Comp Beads are supplied in buffer containing stabilizer and 0.05% sodium azide.
Storage	Store protected from light at 2–8 °C. Do not freeze. The expiration date is indicated on the vial label.

1.1 Background information

Most small molecule and protein-based fluorochromes have broad emission spectra which partially overlap with the spectra of other fluorochromes in flow cytometric analysis. This overlap should be detected and corrected for proper analysis.

The MACS Comp Bead Kit, anti-REA has been developed for optimal compensation of fluorescence spillover of fluorochrome-conjugated REAfinity™ Antibodies. After staining with fluorochrome-conjugated REAfinity Antibodies, the MACS Comp Beads – anti-REA can be used for automated or manual compensation along with the MACS Comp Beads – blank for the control of the negative population.

1.2 Applications

- Compensation of the fluorescence spillover from fluorochrome-conjugated REAfinity Antibodies.

1.3 Reagent and instrument requirements

- MACSQuant® Analyzer, MACSQuant Analyzer 10 (# 130-096-343), or MACSQuant VYB (# 130-096-116)
- MACSQuant Running Buffer (# 130-092-747)
- MACS MiniSampler Plus (# 130-105-745)
- Chill 5 Rack (# 130-092-951)
- MACSQuant Calibration Beads (# 130-093-607)
- Fluorochrome-conjugated REAfinity Antibodies. For more information about antibodies refer to www.miltenyibiotec.com/antibodies.

2. Protocol

2.1 Immunofluorescent staining

1. Mix the MACS Comp Beads – anti-REA and MACS Comp Beads – blank well before use.
2. For each fluorochrome-conjugated REAfinity Antibody used in the experiment label a separate 5 mL sample tube.
3. Add 100 µL of the MACSQuant Running Buffer to each sample tube.
4. Add the appropriate amount of fluorochrome-conjugated REAfinity Antibodies to the appropriate sample tube:
Antibody dilution 1:10 or 1:11: **20 µL**.
Antibody dilution 1:50: **4 µL**.
▲ **Note:** For other formats use 1 µg/mL.
5. Add one full drop of the MACS Comp Beads – anti-REA and one full drop of the MACS Comp Beads – blank to each tube.
▲ **Note:** One full drop is approximately 50 µL.
6. Mix well and incubate for 5–10 minutes in the dark at room temperature. Shake sample during incubation.
7. Dilute each sample by adding 1 mL of the MACSQuant Running Buffer and mix well.
8. (Optional) If no fluorochrome for channel B3 on the MACSQuant Analyzer or B2 on the MACSQuant VYB is used, add 500 µL of MACSQuant Running Buffer and one drop of the MACS Comp Beads – blank to a separate tube and mix well. This sample will be used to set the compensation values for propidium iodide (PI).
▲ **Note:** Fluorochromes for channel B3 or B2 are, for example, PerCP, PE-Cy* 5, or PerCP-Vio* 700.
9. Proceed with calibration of the MACSQuant Instrument (2.2).

2.2 Calibration of the MACSQuant Instrument

▲ Refer to the respective user manual for instructions on how to use the MACSQuant Instrument.

1. Login and prime the MACSQuant Instrument. Check that the system status bar indicates 'Ready' (green).
2. Ensure that the single tube holder is correctly attached to the instrument.
3. Resuspend the MACSQuant Calibration Beads by vortexing.
4. Click on the Barcode reader icon. Position the barcode label of the MACSQuant Calibration Beads vial in front of the barcode reader.

▲ **Note:** A dialog box will appear that instructs the user to add one drop of calibration beads to a single tube.

5. Place the tube containing one drop of Calibration Beads into the single tube holder and click on **run**. The instrument will then proceed to automatically dilute and mix the bead suspension before performing calibration.

▲ **Note:** The system status bar will indicate 'processing sample' (blue) until completion of the task.

6. Following successful calibration the status bar will report 'MACSQuant ready; Calibration ok'; this is visually represented by a green status.

▲ **Note:** The gain and trigger settings are now preset for running experiments and can be viewed under the channel tab in the **Instrument Settings**.




7. Proceed to compensation (2.3).

2.3 Compensation

1. Change forward scatter (FSC) and side scatter (SSC) instrument settings from lin to **log3** and reduce the SSC voltage by 50 V.
2. Click on the **Experiment** tab and choose **Chill 5 tube rack**. Place the sample tubes into the rack.
3. Select the chosen rack positions and group them by clicking **Group** in the rack template.
4. Choose **Express** mode in the **Settings** tab. Choose type **Setup** and mode **CompensationMultiColor** from the dropdown menus.
5. Choose the appropriate color for each rack position under **SampleID**.
6. (Optional) If using MACS Comp Beads – blank for PI compensation, select **PI** in the **SampleID** for this tube.
7. Click on run and follow the given instructions by the MACSQuant Instrument.

3. Example of a compensation using the MACS Comp Bead Kit

The MACS Comp Beads – anti-REA were stained with CD52 antibodies conjugated to FITC and PE (clone REA164) and analyzed by flow cytometry using the MACSQuant Analyzer. Gating was performed on single bead events (A). MACS Comp Beads are shown before (B) and after (C) compensation using the Compensation Multicolor program with the MACSQuant Analyzer.

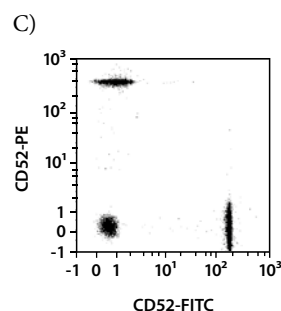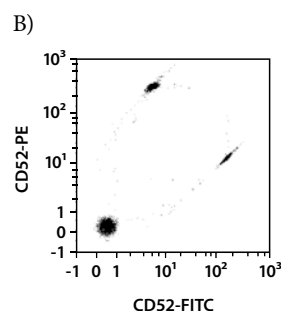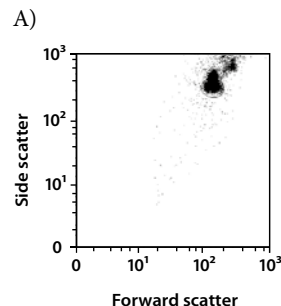

Refer to www.miltenyibiotec.com for all data sheets and protocols. Miltenyi Biotec provides technical support worldwide. Visit www.miltenyibiotec.com for local Miltenyi Biotec Technical Support contact information.

Legal notices

Limited product warranty

Miltenyi Biotec B.V. & Co. KG and/or its affiliate(s) warrant this product to be free from material defects in workmanship and materials and to conform substantially with Miltenyi Biotec's published specifications for the product at the time of order, under normal use and conditions in accordance with its applicable documentation, for a period beginning on the date of delivery of the product by Miltenyi Biotec or its authorized distributor and ending on the expiration date of the product's applicable shelf life stated on the product label, packaging or documentation (as applicable) or, in the absence thereof, ONE (1) YEAR from date of delivery ("Product Warranty"). Miltenyi Biotec's Product Warranty is provided subject to the warranty terms as set forth in Miltenyi Biotec's General Terms and Conditions for the Sale of Products and Services available on Miltenyi Biotec's website at www.miltenyibiotec.com, as in effect at the time of order ("Product Warranty"). Additional terms may apply. BY

USE OF THIS PRODUCT, THE CUSTOMER AGREES TO BE BOUND BY THESE TERMS.

THE CUSTOMER IS SOLELY RESPONSIBLE FOR DETERMINING IF A PRODUCT IS SUITABLE FOR CUSTOMER'S PARTICULAR PURPOSE AND APPLICATION METHODS.

Technical information

The technical information, data, protocols, and other statements provided by Miltenyi Biotec in this document are based on information, tests, or experience which Miltenyi Biotec believes to be reliable, but the accuracy or completeness of such information is not guaranteed. Such technical information and data are intended for persons with knowledge and technical skills sufficient to assess and apply their own informed judgment to the information. Miltenyi Biotec shall not be liable for any technical or editorial errors or omissions contained herein.

All information and specifications are subject to change without prior notice. Please contact Miltenyi Biotec Technical Support or visit www.miltenyibiotec.com for the most up-to-date information on Miltenyi Biotec products.

Licenses

This product and/or its use may be covered by one or more pending or issued patents and/or may have certain limitations. Certain uses may be excluded by separate terms and conditions. Please contact your local Miltenyi Biotec representative or visit Miltenyi Biotec's website at www.miltenyibiotec.com for more information.

The purchase of this product conveys to the customer the non-transferable right to use the purchased amount of the product in research conducted by the customer (whether the customer is an academic or for-profit entity). This product may not be further sold. Additional terms and conditions (including the terms of a Limited Use Label License) may apply.

CUSTOMER'S USE OF THIS PRODUCT MAY REQUIRE ADDITIONAL LICENSES DEPENDING ON THE SPECIFIC APPLICATION. THE CUSTOMER IS SOLELY RESPONSIBLE FOR DETERMINING FOR ITSELF WHETHER IT HAS ALL APPROPRIATE LICENSES IN PLACE. Miltenyi Biotec provides no warranty that customer's use of this product does not and will not infringe intellectual property rights owned by a third party. BY USE OF THIS PRODUCT, THE CUSTOMER AGREES TO BE BOUND BY THESE TERMS.

Trademarks

MACS, MACSQuant, the Miltenyi Biotec logo, REAfinity, and Vio are registered trademarks or trademarks of Miltenyi Biotec and/or its affiliates in various countries worldwide.

Cy is a registered trademark of GE Healthcare UK Limited.

Copyright © 2021 Miltenyi Biotec and/or its affiliates. All rights reserved.